


# ALLEVYN Gentle Border

## Advanced Foam Wound Dressings

Proud to Put People First  
**ALLEVYN**◊  
Classic

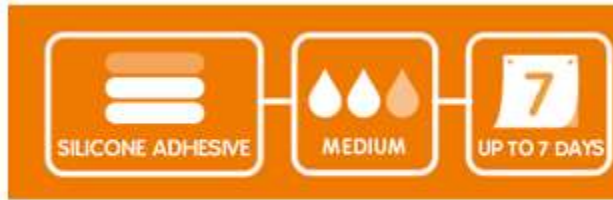

### ALLEVYN Gentle Border

ALLEVYN Gentle Border Range includes:

- ALLEVYN Gentle Border
- ALLEVYN Gentle Border Multisite
- ALLEVYN Gentle Border Sacrum
- ALLEVYN Gentle Border Heel



The ALLEVYN Gentle Border range has been designed for people with particularly sensitive or fragile skin.<sup>1,2,3</sup> These dressings have a soft silicone gel adhesive which minimizes trauma to the wound at dressing change,<sup>1,3,4</sup> helping to avoid patient pain without compromising care.<sup>1,3</sup>

- ALLEVYN triple-action technology maintains the optimal fluid balance to maintain a moist wound environment.<sup>5</sup>
- Conformable dressing – providing comfort for the patient.<sup>1,3</sup>
- Minimizes trauma to the wound and pain to the patient during dressing change.<sup>1,3</sup>
- Easy to apply and remove<sup>1,3,6</sup> – the dressing can be repositioned easily without losing its adherent properties.<sup>1</sup>
- Showerproof dressing – allowing patients to shower with the dressing in place.<sup>7</sup>
- Available in a range of shapes and sizes – squares, rectangles and a range of shapes for dressing awkward areas.<sup>7</sup>

We have collaborated with clinicians and patients to develop a unique range of ALLEVYN Gentle Border Dressings designed especially for the contours of the human body.

ALLEVYN Gentle Border Dressings are intended for use on:

- Shallow, granulating wounds
- Chronic and acute exudative wounds

- Full- and partial-thickness wounds such as pressure ulcers, leg ulcers and diabetic foot ulcers
- Infected wounds
- Malignant wounds
- Surgical wounds
- First and second degree burns
- Donor sites




#### ALLEVYN Gentle Border Size and Ordering Information

Product Name	Product Number	Size	Pieces per Carton	Cartons per Case	HCPCS Code
ALLEVYN Gentle Border	66800276	3" x 3"/7.5cm x 7.5cm	10	4	A6212
ALLEVYN Gentle Border	66800270	4" x 4"/10cm x 10cm	10	10	A6212
ALLEVYN Gentle Border	66800279	5" x 5"/12.5cm x 12.5cm	10	4	A6212
ALLEVYN Gentle Border	66800280	7" x 7"/17.5cm x 17.5cm	10	4	A6213
ALLEVYN Gentle Border Multisite	66800959	6 3/4" x 7 1/16"/17.1cm x 17.9cm	10	6	A6212
ALLEVYN Gentle Border Sacrum	66800898	6 5/8" x 6 3/4"/16.8cm x 17.1cm	10	6	A6212
ALLEVYN Gentle Border Sacrum	66801031	8 1/2" x 9"/21.6cm x 23cm	10	8	

ALLEVYN Gentle Border Heel	66800506	9" x 9 1/8"/23cm x 23.2cm	5	7	A213
----------------------------	----------	---------------------------	---	---	------

Like all ALLEVYN Dressings, ALLEVYN Gentle Border has at its heart a unique triple layer technology that effectively manages exudate, helping to ensure that the wound is kept at optimum moisture levels to support moist wound healing.<sup>5</sup>

The dressing consists of a highly absorbent<sup>8</sup> hydrocellular foam pad held between a perforated wound contact layer which is coated with a soft silicone gel adhesive and a highly permeable outer top film.<sup>9,10,11</sup>

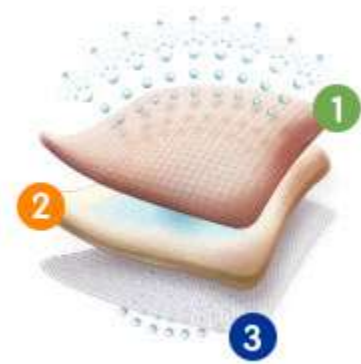


ALLEVYN Gentle Border Dressings achieve the benefits of easy application and removal,<sup>1,3,6</sup> and the ability to reposition through an adhesive technology based on a soft silicone gel.

ALLEVYN Gentle Border Dressings, with their silicone gel adhesive, are designed to minimize trauma to the wound and pain to the patient at dressing changes.<sup>1,3</sup>


In a volunteer evaluation, pain on dressing removal was recorded at day 1 and day 3.<sup>1</sup>

Pain on Removal of ALLEVYN Gentle Border and Competitor Dressing, 1 and 3 Days Wear

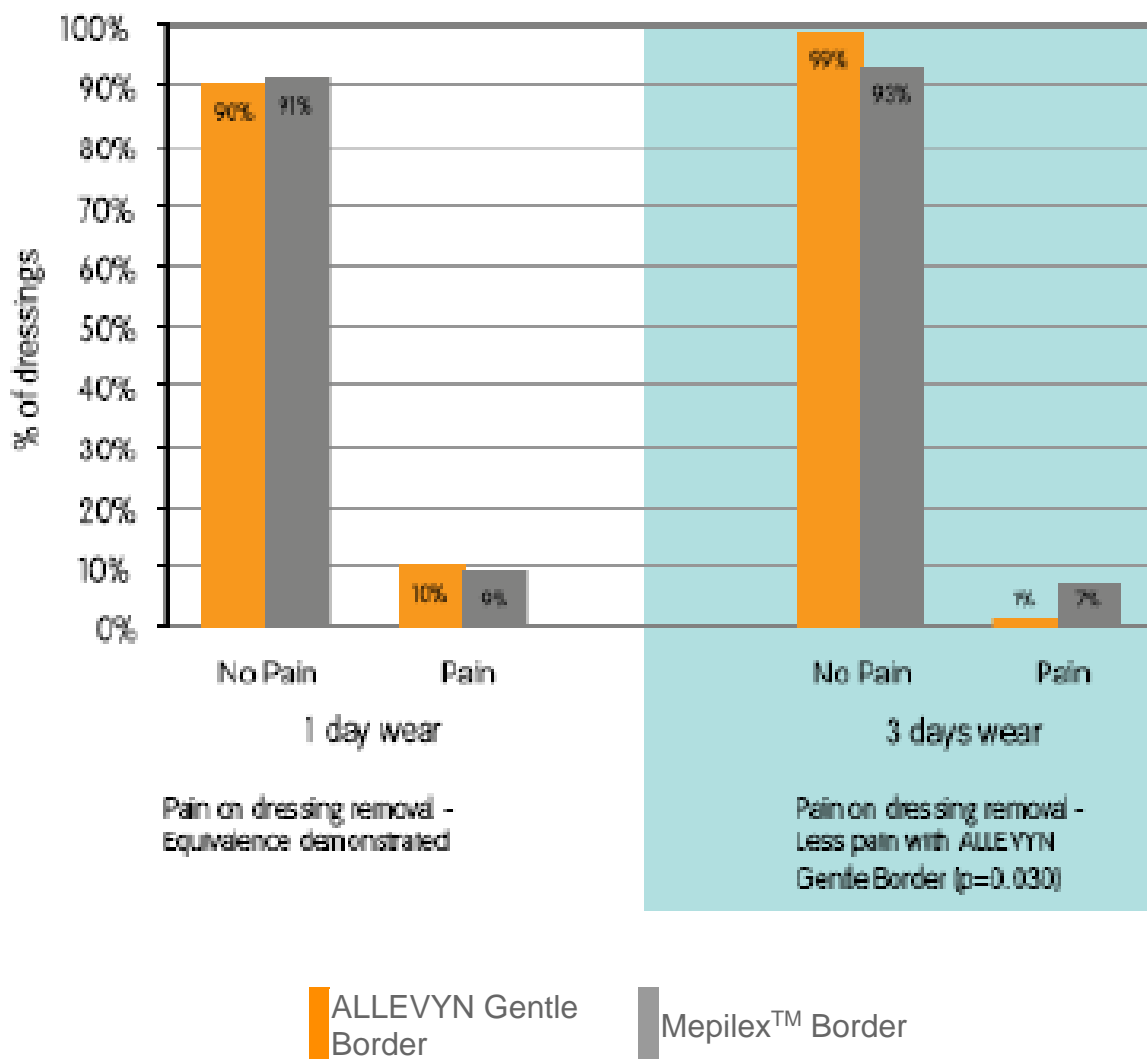

In a volunteer evaluation more ALLEVYN Gentle Border Dressings were still in place after 1,3 and 7 days compared to Mepilex Border<sup>1</sup>

Percentage of Dressings in Place of ALLEVYN Gentle Border and Competitor Dressings after 1,3 and 7 Days Wear

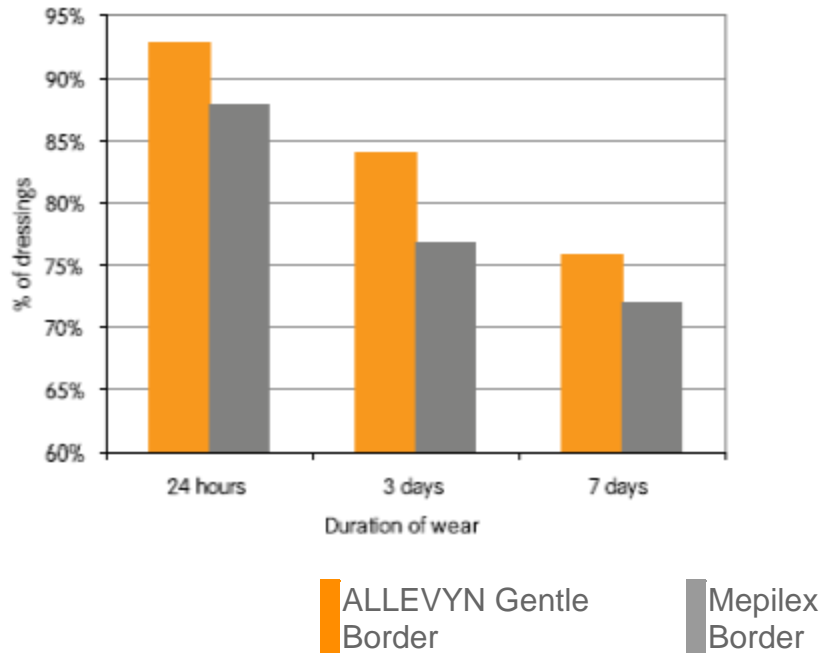

References:

1. S. Brown, S. Palmer. An open, prospective, randomized, within volunteer comparison of ALLEVYN Gentle Border and Mepilex Border. Data on file: OR-DOF/001. October 2008
2. S. W. Ohen. M.D. "Ageing changes in skin," <http://nim.nih.gov/medlineplus/ency/article/004014.htm> update on 22nd August 2006
3. T. Hurd et al. A multicentre in-market evaluation of ALLEVYN Gentle Border. Wounds UK, 2009, Vol 5. No 3
4. G. Walton. ALLEVYN Gentle Border, A 17 day wound healing study in full thickness wounds in minipigs. October 2008
5. D. Myers. The clinical and physical properties of ALLEVYN Gentle Border dressing and the moist wound environment. November 2007
6. S. Palmer, C. Mistry. An open, prospective, randomised, comparative study to compare the performance of ALLEVYN Gentle Border Multisite with an alternative silicone adhesive dressing. Data on file report: OR-DOF/017. June 2011
7. D. Hussein. Results from a multi-centre, non-comparative clinical in market evaluation of ALLEVYN Gentle Border dressings (CIME/008) and from open, prospective, randomised within volunteer (VTE/HVT 050) comparison of Dakota and Mepilex Border in terms of showerproof and use on awkward areas. Data on file: OR-DOF/011. August 2010
8. L. Tompkins. Data on File Report. 0711080. November 2007
9. N. Patrick. Technical Justification Summary AGB007. October 2010

10. S. Carpenter, D Foster. Data on File Report: DS/07/187/R1v2. November 2007

11. S. Carpenter, D Foster. Data on file Report: DS/08/043/Ra. June 2008