

Progressive Resting Splints

Keep the gains you've made!

Turnbuckles for every extremity



Ankle, p. 311



Knee, p. 341



Elbow, p. 377



Turnbuckle Functional Position Splints

These static progressive splints help maintain gains in ROM achieved through stretching and splint wear. The steel rod-and-turnbuckle mechanism provides a low-load stretch to the joint to gradually reduce a wrist flexion contracture while maintaining a functional hand position. Turnbuckle with triangular-shaped grip adjusts easily with the splint in place for optimal positioning.

Rigid Kydex® splint shell is hinged at the wrist, with the turnbuckle connecting the hand and forearm portions of the splint. Soft padded two-part Beta Pile II™ cover is comfortable against skin and is removable for washing. Padded straps are reinforced with hook-and-loop closures in areas of repetitive closure for added durability.

Splint adjusts from 80° of wrist flexion to 60° of wrist extension. For appropriate fit, patient must have passive finger extension, and wrist contracture cannot exceed 90° flexion.

Kydex splint shell • Beta Pile™ liner: removable and washable
• Heat-adjustable

Sizing is approximate. B fits smaller Wmn; C fits larger Wmn and smaller Men; D fits most Men.



Triangular-shaped grip for easier adjustment

SIZE	LEFT	RIGHT	QTY
B	#2970002740	#2970002743	EACH
C	#2970002741	#2970002744	EACH
D	#2970002742	#2970002745	EACH

Must be used under a physician's supervision.

IMP • CA California residents, see left page.


Suggested code: L3931

Adjust each finger individually...

Shown in Light Blue



Shown in Dark Blue

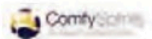

Hand/Wrist Separate Finger Orthoses

support wrist and whole hand, allowing custom adjustment of each finger individually. Finger separator included. Adjustable metal frame bends to fit without heat or additional tools. Slim, nonremovable, surface-cleanable Comfyprene™ neoprene cover. One size.

Sizing: Measure wrist to tip of middle finger. 6"-8"; 12½"L. Pediatric sizes on p. 403.

COLOR	LEFT	RIGHT	QTY
DARK BLUE	#2970002447	#2970002505	EACH
LIGHT BLUE	#2970002448	#2970002503	EACH

CA California residents, see left page.

Flexible wrist and fingers and progressive thumb abduction

DynaPro™ Finger Flex Orthoses

Flexible wrist, MCP, and PIP joints with flexible finger platform designed for treatment of spasticity. "Flex" properties will "give" during involuntary muscle contractions to help inhibit tone.

Resting Hand Orthoses (Thumb Ease) combine the overall design of the Finger Flex (including the flexible wrist segment) with the rigid finger platform of a traditional resting hand splint to help manage wrist spasticity while maintaining functional hand alignment.

Both models have a thumb pocket with two gel inserts to help progressively increase thumb abduction. Heat-moldable core for custom fitting. Finger separator included. Color-coded straps for easy application and staff training. Antibacterial cover helps reduce odor. Includes extra cover for laundering.

Sizing: Measure width of MP joint and distance from wrist to tip of middle finger. See chart below.

SIZE	WIDTH OF MP JOINTS	WRIST-TO-FINGERTIP LENGTH
SMALL	UP TO 3"	UP TO 6"
MEDIUM	3"-3½"	6½"-8"
LARGE	3½"-4"	8"-9½"



Color-coded straps



Finger separator included

SIZE	LEFT	RIGHT	QTY
FINGER FLEX			
SMALL	#2970002700	#2970002699	EACH
MEDIUM	#2970002702	#2970002701	EACH
LARGE	#2970002704	#2970002703	EACH
RESTING HAND (THUMB EASE)			
SMALL	#2970002706	#2970002705	EACH
MEDIUM	#2970002708	#2970002707	EACH
LARGE	#2970002710	#2970002709	EACH

IMP • CA California residents, see left page.

LOG IN FOR PRICING at AliMed.com

UPPER EXTREMITIES

PROGRESSIVE RESTING SPLINTS

FAX 800.437.2966

866.662.0283 AliMed.com

369